

PHILIPS

IPL hair removal system

Lumea Comfort

For use on body

22 minutes to treat lower legs

Lifetime >100,000 light pulses

Compact corded design

SC1981/50



Prevents hair regrowth. Provides up to 75% hair reduction*

*when used on legs bi-weekly

Philips Lumea Comfort IPL hair removal system works wonders to prevent the reappearance of hair on your body. Gentle pulses of light, applied regularly, keep skin silky-smooth every day.

Long lasting smoothness

- Effortlessly effective
- Enjoy smooth skin every day

Hair removal with IPL technology

- For effective use in the safety and comfort of your home

Gentle and effective on body


- Five adjustable light energy settings
- Gentle and effective even on sensitive areas
- Skin Color Sensor
- Slide and Flash mode for a more convenient application
- Suitable for different hair and skin types

Most compact Lumea design

- No replacement parts, no hidden costs


Highlights

Enjoy smooth skin every day

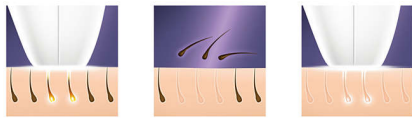

Philips Lumea applies gentle pulses of light to the hair root. This results in the hair shedding naturally and hair regrowth is inhibited. Repeating this treatment regularly will keep your skin smooth everyday.

Effortlessly effective

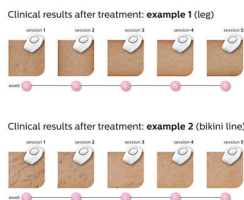

Our clinical studies have shown significant hair reduction in as little as four bi-weekly treatments resulting with smooth hair free skin. To maintain these results, simply repeat the treatment when needed. Time between treatments may vary based on your individual hair regrowth. - To enlarge this picture, please click on the image in the image gallery on the top of this page

IPL technology

Philips Lumea Comfort uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in professional beauty salons. Philips has adapted this technology for effective use in the comfort of your own home. Philips worked closely with skin care experts to develop this technology. For

more than 10 years Philips carried out extensive consumer research on IPL technology with more than 2000 women.

Suitable hair and skin types

Haarfarbe Haarkleur Hair color Couleur de poils	Hautton / Huidskleur / Skin tone / Couleur de la peau	1	2	3	4	5	6
Dark brown / Donkerbruin / Dark brown / Brun foncé		✓	✓	✓	✓	✗	✗
Medium brown / Lichtbruin / Medium brown / Brun moyen		✓	✓	✓	✓	✗	✗
Light brown / Lichtblond / Light brown / Blond clair		✓	✓	✓	✓	✗	✗
Dark blonde / Donkerblond / Dark blonde / Blond foncé		✓	✓	✓	✓	✗	✗
Light blonde / Lichtblond / Light blonde / Blond clair		✗	✗	✗	✗	✗	✗
Red / Rood / Red / Rouge		✗	✗	✗	✗	✗	✗
White / Wit / White / Blanc		✗	✗	✗	✗	✗	✗
Grey / Grijs / Grey / Gris		✗	✗	✗	✗	✗	✗

Philips Lumea Comfort IPL Hair Removal System works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Lumea Comfort IPL Hair Removal System is NOT effective on red, light-blonde or white/grey hair. - To enlarge this picture, please click on the image in the image gallery on the top of this page

Application areas

Use with confidence to remove hair on the body including legs, arm pits, bikini area, belly and arms.

5 light energy settings



Philips Lumea IPL Hair Removal System has five adjustable light energy settings to ensure gentle but effective treatment. Used correctly, Philips


Lumea's IPL treatment is comfortable and gentle to use even on sensitive skin and sensitive body areas.

Skin Color Sensor



The unique Skin Color Sensor measure the treated skin complexion at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for Lumea Comfort, it will automatically stop emitting pulses.

Slide and Flash mode

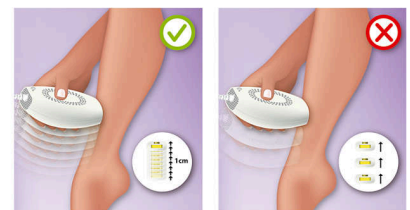

Lumea Comfort IPL Hair Removal System enables a sliding motion for a more convenient and continuous application.

No replacement parts needed

Philips Lumea Comfort IPL Hair Removal System is a complete solution that comes ready to work straight out of the box. More importantly, it does not require any replacement lamps, providing 4 years of full body treatments with no extra cost

Specifications

Application areas

Body areas: Legs, Belly, Bikini, Armpits

Sensitive areas: small areas such as armpit

Safety and adjustable settings

Integrated UV filter: Protects skin from UV light

5 light energy settings: Adjustable to your skin type

Integrated safety system: Prevent unintentional flashing

Skin Tone sensor: Detects your skin tone

Application mode

Slide and Flash: For easy motion application

Step and Flash: For treatment on curvy areas

Corded / cordless use: Corded

Power

Voltage: 100-240 V

Technical specifications

High performance lamp: Delivers <gt/> 100,000 flashes

Service

Warranty: 2 years global warranty

Items included

Instructions for use: Quick Start Guide, User manual

Adapter: 24V / 1500mA

Application time

Lower leg: 11 min

Armpit: 1,5 min

Bikini line: 1,5 min

Sensitive areas: 2,5 min

